



BICRON ♦ NE
6801 Cochran Road
Solon, Ohio 44139
Telephone: (216) 248-7400
(800) 472-5656
Telefax: (216) 349-6581

Surveyor 2000™
Portable Survey Meter

User's Manual

Publication No. 1016-0-U-0896-001

* * * Release Date * * *

August 30, 1996

Part No. 1016907
Rev. G

Division of Saint-Gobain/Norton Industrial Ceramics Corporation

ISO 9001 Quality Certified



Surveyor 2000™
Portable Survey Meter

User's Manual

Publication No. 1016-0-U-0896-001

*** * * Release Date * * ***

August 30, 1996

Part No. 1016907
Rev. G

FOREWORD

This manual provides the basic operation and maintenance procedures for the BICRON Surveyor 2000™ Portable Survey Meter.

Section 1.0 Introduction provides a general description of the instrument and its operation, and a detailed listing of its physical and performance specifications.

Section 2.0 Battery Installation describes the procedure for changing the battery and checking its performance.

Section 3.0 High Voltage Test describes how to read the High Voltage.

Section 4.0 Operation provides a detailed description of all the features and settings for the Surveyor 2000, and complete operating instructions.

Section 5.0 Hardware Description provides a brief description of the three electronic circuits that make up the instrument, and a description of the GM Tube.

Section 6.0 Calibration provides directions for calibration of the instrument.

The **Appendices** are: A) QC Acceptance Procedure which includes calibration procedures, B) a complete spare parts list so instruments can be repaired on-site, and C) schematic and pictorial diagrams to facilitate repair procedures.

Writing Conventions

In order to maintain consistency throughout this and all BICRON ♦ NE manuals, certain writing conventions have been followed for safety warnings. They are divided into

three categories and defined as follows:

- **DANGER** indicates an imminently hazardous situation which, if not avoided, will result in death or serious injury. ***DANGER NOTICES ALWAYS APPEAR IN BOLD, ITALICIZED UPPER CASE LETTERS.***
- **WARNING** indicates a potentially hazardous situation which, if not avoided, could result in death or serious injury. ***WARNING NOTICES ALWAYS APPEAR IN UPPERCASE BOLD LETTERS.***
- **CAUTION** indicates a potentially hazardous situation which, if not avoided, may result in minor or moderate injury. It may also be used to alert against unsafe practices. ***CAUTION notices always appear in bold, italicized letters.***

The definition of these safety warnings is according to ANSI Z535.4. The style of the warnings (bold, italicized, etc) is BICRON ♦ NE's.

In addition to the above, we have added the following warning:

- **NOTE** indicates a situation which has the potential for erroneous data collection, loss of electronic data, or damage to equipment, but which does not directly affect the safety of the operator with respect to this product. The responsibility for any safety consequences as a result of erroneous data lies solely with the operator. ***NOTE notices always appear in italics.***

TABLE OF CONTENTS

Section	Page
BICRON ♦ NE Warranty Statement	5
Procedures and Cautions	6
1.0 Introduction	
1.1 General Description	7
1.2 Specifications	7
2.0 Battery Installation	9
2.1 Procedure	9
2.2 Test	9
3.0 High Voltage Test	9
4.0 Operation	
4.1 Radiation Measurements	10
4.2 Selectable Response Time	10
4.3 Audio	10
5.0 Hardware Description	
5.1 Circuit Description	11
5.2 GM Tube (Probes)	11
6.0 Calibration	12
6.1 Standard (mR/H)	12
6.2 Optional (CPM) Calibration	12
Appendices	
A QC Acceptance Procedure, Number 1016937	13
B Spare Parts List, Number 1016917	15
C Drawings	19

**BICRON ♦ NE WARRANTY STATEMENT
COVERING PORTABLE MONITORS**

Instruments and options manufactured by BICRON ♦ NE are warranted against defects in materials and workmanship for a period of two years from the date of shipment, unless otherwise agreed upon by BICRON ♦ NE and the customer in writing.

BICRON ♦ NE's obligation with regard to such products shall be limited to repair or replacement FOB BICRON ♦ NE factory or authorized repair station, at BICRON ♦ NE's option.

The calibration (when applicable) for each system is warranted to be within its specified accuracy at the time of shipment. If this initial calibration is determined to be in error, the system will be recalibrated at no charge.

The aforesaid warranty does not cover systems, options or probes which are subject to excessive physical abuse or are used for purposes other than those intended. In no event shall BICRON ♦ NE be liable for consequential or special damages, transportation, installation, adjustment, work done by customer, or other expenses which may arise in connection with such defective product or parts.

EXCLUSION OF LIMITED WARRANTY

THERE ARE NO WARRANTIES, EXPRESS OR IMPLIED, INCLUDING WITHOUT LIMITATION ANY IMPLIED WARRANTY OF MERCHANTABILITY OR FITNESS, WHICH EXTEND BEYOND THE DESCRIPTION ON THE FACE HEREOF. THIS EXPRESS WARRANTY EXCLUDES COVERAGE OF, AND DOES NOT PROVIDE RELIEF FOR, INCIDENTAL OR CONSEQUENTIAL DAMAGES OF ANY KIND OR NATURE, INCLUDING, BUT NOT LIMITED TO, LOSS OF USE, LOSS OF SALES OR INCONVENIENCE. THE EXCLUSIVE REMEDY OF THE PURCHASER IS LIMITED TO REPAIR, RECALIBRATION, OR REPLACEMENT OF THE SYSTEM AT BICRON ♦ NE's OPTION.

This warranty specifically excludes the following items which are covered by their original manufacturers' warranties: photomultiplier tubes, GM and proportional tubes, crystal and other solid-state detectors and batteries.

PROCEDURES and CAUTIONS

The equipment herein described is designed and manufactured in compliance with all applicable safety standards. Nevertheless, certain hazards are inherent in the use of electronic and radiometric equipment.

Adequate warnings are included in the manual and on the product itself to cover hazards that may be encountered in normal use and servicing of this equipment. No other procedures are warranted by Bicon.

It shall be the owner's or user's responsibility to see to it that the procedures and cautionary notes are heeded.

Failure on the part of the user in any way to follow the prescribed procedures shall absolve Bicon and its agents from any resulting liability.

This instrument is intended solely for the detection and measurement of ionizing radiation. It should be used only by persons who have been trained in the proper interpretation of its readings and the appropriate safety procedures to be followed in the presence of radiation.

All instructions and warnings contained in this manual or on the instrument must be read before use and must be strictly followed. Failure to follow these instructions and warnings may result in inaccurate readings and/or user hazard.

Indicated battery and other operations tests must be performed prior to each use to assure that the instrument is functioning properly.

CAUTION

FAILURE TO CONDUCT PERIODIC PERFORMANCE TESTS IN ACCORDANCE WITH ANSI N323-1978, PARAGRAPHS 4.6 AND 5.4, AND TO KEEP RECORDS THEREOF IN ACCORDANCE WITH PARAGRAPH 4.5 OF THE SAME STANDARD, COULD RESULT IN ERRONEOUS READING OF POTENTIAL DANGER. ANSI N323-1978 BECOMES, BY THIS REFERENCE, A PART OF THIS OPERATING PROCEDURE.

INSPECTION

Instruments should be examined and tested as soon as received. Claims for transportation damages, if any, should be filed at once with the delivery carrier.

1.0 Introduction**1.1 General Description**

The Bicon Surveyor 2000 is a portable survey meter designed for use with an appropriate GM probe to detect and measure ionizing radiation.

The instrument features an internal, energy-compensated GM detector (0-2000 mR/h range only), a recessed meter movement, laminated control panel, cpm and mR/h meter scales, single on-off/range selector switch, MHV probe connector and mounted probe holder.

Advanced circuit design provides a detector HV check, anti-saturation, dead time compensation, switch selectable response time (optimized for each range), and built in audio.

1.2 Specifications**Radiation Detected:**

Alpha (depending on probe used); beta; gamma with external probe; gamma and x-ray with internal detector

Detector:

GM tube, internal; choice of GM probes, external

Range:

Linear ranges of 0-0.2, 0-2, 0-20, 0-200 mR/h (SWGGM or EWGM probe) 0-2000 mR/h (internal detector only); 0-240; 0-2400; 0-24,000; 0-240,000 CPM (external probes)

Accuracy:

Within 10% of reading for ^{137}Cs when calibrated according to NRC Reg. Guide 10.8

High Voltage:

Electronically stabilized, factory set at 900 V

HV Test:

Exclusive self test to verify detector HV power supply

Connector: MHV**Warmup time:** None**Saturation:**

Typically greater than 1000 R/h on all ranges (with exclusive anti-saturation circuit) for internal detector and most GM probes; greater than 5 R/h for pancake GM probes

Response Time:

Switch-selectable, optimized for each range. 0-90% of final reading as follows:

Range	Time	
	Fast	Slow
X0.1	6 sec.	25 sec.
X1	2 sec.	6 sec.
X10	1 sec.	3 sec.
X100	< 1 sec.	1 sec.
X1000	< 1 sec.	1 sec.

Dead time Compensation:

Exclusive circuitry provides a near linear response.

Temperature:

Operational from -40°C to $+60^{\circ}\text{C}$

Humidity:

less than 5% change in reading from 10 - 95% relative humidity

Control:

Eight position rotary switch: "off", "bat", "HV", "X1000", "X100", "X10", "X1", "X0.1"; rotary response time switch

1.0 Introduction (cont'd)**1.2 Specifications (cont'd)****Battery Complement:**

Single 9 Volt (MN1064 or equal). The additional battery holder may be used for storage of spare or parallel-wired

Battery Life:

Greater than 100 hours with single battery or greater than 200 hours with parallel batteries option

Display:

Ruggedized, recessed, high-torque 1 mA meter with 8.51 cm. (3.35 inch) scale marked "0-2 mR/h", "0-2400 CPM", "bat. ok", "HV ok". Meter protected by impact-resistant Lexan® polycarbonate window

Geotropism:

Within $\pm 2\%$ of full scale

Shock:

100 g per lightweight machine of MIL-STD 202C, method 202B

Vibration:

5 g in each of three mutually orthogonal axes at one or more frequencies from 10-33 Hz

Construction:

Splash-proof, shock-proof, two-piece, all-metal case. Scratch-resistant laminated control panel and Bicon Kleen-Krome® trim on case top, durable black polyurethane paint on handle and case bottom, stainless steel probe holder.

Audio:

A built in speaker, controlled by a panel mounted on-off switch, provides an audible "click" for each detector pulse or (if desired) an audible over-range alarm signal.

Size:

10.8 x 20.3 x 17.3 cm (4.25" x 8" x 6.8") including handle and probe clip

Weight: 1 kg. (2.2 lb), excluding probe

2.0 Battery Information

Battery Type: 9-volt Mallory MN1604 or equivalent.

2.1 Installation Procedure

1. Turn instrument off.
2. Open pull catches at ends of case and separate case bottom from top.
3. Install battery in appropriate clip on bottom circuit board (clip for spare battery is so marked), observing proper polarity.
4. Replace bottom, orienting rubber under battery; close catches.

Parallel-wired Option: Instruments with this option have both battery clips wired into the circuit. Installing a second battery thus provides approximately twice the operational hours of one battery. Only one battery is needed to power the instrument, however.

2.2 Battery Test

To check the condition of the battery turn the Control Switch to the "bat." position (Figure 2). If the battery has sufficient charge, the meter will read within the "bat. ok" range (Figure 1).

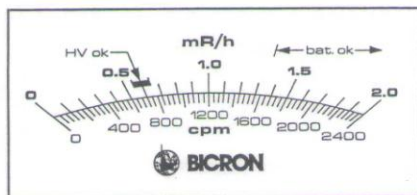


Figure 1
Surveyor 2000 Meter Scale

3.0 High Voltage Test

To check the high voltage setting turn the Control Switch to the "HV" position (Figure 2). If the high voltage is correctly set, the meter will read within the "HV ok" range (Figure 1).

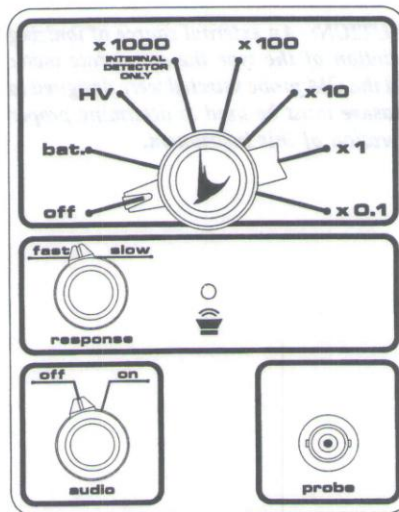


Figure 2
Surveyor 2000 Control Panel

4.0 Operation

4.1 Radiation Measurements

1. When using an external GM probe, select the appropriate response time (Section 4.2 Selectable Response Time) and measurement range (Figure 2). Use only the "X0.1", "X1", "X10", or "X100" ranges. Read only the cpm scale unless using a SWGM (side-wall GM) or EWGM (end window GM) probe which has been calibrated in mR/h with the unit.
2. When using the internal detector, select the appropriate response time and turn the control switch to the X1000 range. Read the mR/h scale only. The detector is mounted near the front of the case, centered vertically and slightly to the right of center horizontally. Position the instrument so that radiation is incident on the front of case.

CAUTION: An external source of ionizing radiation of the type the count rate meter and the GM probe selected were designed to measure must be used to determine proper operation of this instrument.



Figure 1
Surveyor 2000 Count Rate Meter

4.2 Selectable Response Time

In addition to the built in circuitry which optimizes response time for each range, response time can be controlled by the "response" Switch. "slow" or "fast" response time can be selected.

4.3 Audio

An internal speaker produces an audible "click" for each detector pulse. This speaker is controlled by a panel mounted on-off switch labeled "audio".

An audible alarm sounds when the meter is approximately 30% or more above full scale on any range. This over-range alarm can be defeated by turning off an internal on-off switch located on the main circuit board.

The over-range alarm functions only when the panel mounted Audio Control Switch is in the "off" position and the internal switch is "on". When the speaker control is in the "on" position, the over-range alarm is defeated.



Figure 2
Surveyor 2000 Internal Circuit Board

5.0 Hardware Description

5.1 Circuit Description

The electronic circuitry in the Bicon surveyor 2000 is contained on three interconnected printed circuit boards. Modern solid-state integrated circuitry is used throughout. The major components are:

1. The high voltage power supply which is a feedback-regulated, electronically stabilized supply for the GM tube potential. Additional circuitry provides an HV test readout on the meter scale.
2. The count-rate meter which is a linear charge pump rate meter that converts the GM tube pulses to an exposure rate reading on the calibrated meter scale. The circuitry includes a unique dead-time compensation technique to provide nearly linear response over the full range, an anti-saturation circuit which forces the meter beyond full scale in high radiation fields, automatic time constant selection, and temperature compensation.
3. The audio circuitry for individual pulse counting and over-range alarm.

5.2 GM Tube (Probes)

The GM tube consists of a thin shell which acts as the cathode, a fine wire anode suspended within the shell, and an inert gas into which a small amount of a halogen gas has been mixed to act as a quenching agent. End window and pancake tubes have a thin mica entrance window.

A potential of approximately 900 volts is maintained between the two electrodes with the anode always positive. This voltage is slightly less than that required to produce a discharge in the gas. When a nuclear particle or ray of sufficient energy enters the GM tube, it ionizes a molecule of the inert gas. The positive ions are attracted to the cathode and the electrons are attracted to the anode because of the high voltage maintained between the electrodes. In their movement toward the electrodes, these charged particles trigger the ionization of additional gas molecules, resulting in an avalanche of ions flowing between the electrodes. The gas discharge thus created is similar to the glow of a neon lamp. The tube conducts as long as the gas is in the ionized state.

The small amount of halogen gas in the gas mixture quenches the flow of ions, suppressing further electron avalanches until another nuclear particle or ray enters the tube. This flowing and quenching results in a rapid pulse or surge of current in the external circuit. The number of pulses per minute is approximately proportional to the radiation exposure rate. The meter, suitably connected to the tube, indicates the exposure rate or counts per minute (depending on probe used) on a calibrated scale.

6.0 Calibration

This section provides a general description of calibration conditions. Detailed calibration procedures are part of the Q.C. Acceptance Procedure in Appendix A.

6.1 Standard (mR/H)

The instrument is normally factory calibrated using ^{137}C 's gamma rays for the internal detector and when furnished with a specific model SWGM (side wall GM) or EWGM (end window GM) probe. Recalibration is required (a) after servicing, (b) if the GM probe is replaced, and at regular intervals specified by the appropriate regulatory agencies.

The SWGM (with shield closed) or EWGM probe is placed in a known radiation field with the axis of the probe perpendicular to the beam. The front of the case is placed perpendicular to the beam when calibrating the internal detector.

Individual calibration controls are provided for each range, and are used to adjust the meter reading to correspond to the known exposure rate. The locations of these controls are indicated on the main circuit board.

External probes are calibrated on the "X0.1", "X1", "X10", and "X100" scales only. The internal detector is calibrated on the "X1000" scale only.

NOTE: Do not disturb the settings of any controls except those marked "X0.1", "X1", "X10", and "X100" ("X1000" with internal detector only).

Calibration procedures should follow those specified by the appropriate regulatory agencies.

6.2 Optional (CPM) Calibration

Electronic calibration (counts per minute) is the only type that can be done when external probes other than a specific SWGM or EWGM are supplied, or when the instrument is not supplied with a given probe.

When an electronic calibration is performed, the unit is connected to a variable frequency pulse generator. The generator is then set at the frequencies needed to produce 65% and 35% of full scale meter readings (cpm) for each range. Calibration controls are set for 65% of full scale readings, and linearity is checked by going to the 35% readings.



BICRON ♦ NE
6801 Cochran Road
Solon, Ohio 44139
Telephone: (216) 248-7400
(800) 472-5656
Telefax: (216) 349-6581

Surveyor 2000™
Portable Survey Meter

QC Acceptance Procedure

Publication No. 1016-0-Q-0896-001

* * * Release Date * * *

August 30, 1996

Part No. 1016937
Rev. A

Publication Number: 1016-0-Q-0896-001
Part No. 1016937, Rev. A

NOTICE

BICRON ♦ NE reserves the right to make changes to any product to improve reliability, function or design and does not assume any liability arising from either application or use of any product or circuit described herein.

Information presented in this publication has been carefully checked and is believed to be entirely reliable; however, no responsibility is assumed for inaccuracies.

Furthermore, the information does not convey to the user any license under the patent rights of BICRON ♦ NE or others.

Specifications subject to change without notice.

LIMITED RIGHTS LEGEND

The data in this manual is considered "proprietary data." It may be duplicated and used by or on behalf of the customer with the express limitations the "proprietary data" may not be disclosed outside the customer or be used for purposes of manufacture without prior permission of the seller.

This legend shall be marked on any reproduction of this data in whole or in part.

© Copyright Saint-Gobain/Norton Industrial Ceramics Corporation, 1996
All rights reserved

Printed in U.S.A.

1016-0-Q-0896

QC Acceptance Procedure

1.0 Preparation

1. Perform a visual inspection of finished product.
2. Remove all 9V batteries and connect a $9.30\text{ V} \pm 0.05\text{ V}$ power source across the main battery terminals on the battery board.

2.0 Calibration

Perform the following calibrations:

1. Turn the Control Switch to "OFF". Mechanically zero the meter via the rear zero adjustment screw on the meter barrel.
2. Turn the Control Switch to "X1". Check the +5 supply at pin 1 of U6 (ICL 7663). The reading should be $5\text{ VDC} \pm 10\%$.
3. Leave the Control Switch set at "X1" and connect a voltmeter between pins 1 (ground) and 15 of the 24 pin connector. Adjust R31 (50Kohm zero pot) so the voltmeter reads $1\text{ mVDC} +4/-1$.
4. Turn the Control Switch to "HV" and

connect a high voltage measuring device with an impedance ≥ 1000 megohms to the high voltage connection on the internal GM tube. Adjust the high voltage supply via R52 (50 Kohm pot) to $+900\text{ VDC} \pm 3\%$.

5. Connect a high voltage measuring device to the external probe connector and rotate the Control Switch through all 8 positions. You should detect HV only at the probe connector in "X100", "X10", "X1" and "X0.1" positions.
6. Connect the HV measuring device across the internal GM tube and rotate the Control Switch through all 8 positions. You should detect HV in "BAT.", "HV", and "X1000" only.
7. Turn the Control Switch to HV and adjust R34 (500 ohm span pot) until the meter reads in the center of the HV OK checkband.
8. Perform a CPM calibration (unless the instrument is to be supplied with a specific SWGM or EWGM probe):

Table T-1

Range	Pulse Generator output (Hz)	Simulated CPM	Acceptable Meter Reading
X100 (65%)	1643	156,000	140,400 - 171,600
X100 (35%)	1066	84,000	75,600 - 92,400
X10 (65%)	246	15,600	14,040 - 17,160
X10 (35%)	136	8,400	7,560 - 9,240
X1 (65%)	26	1,560	1,404 - 1,716
X1 (35%)	14	840	756 - 924
X0.1 (65%)	2.6	156	140 - 171.6
X0.1 (35%)	1.4	84	75.6 - 92.4

NOTE: The values in this table may not correspond to those for other Bicon or competitive models.

QC Acceptance Procedure

2.0 Calibration (cont'd)

- reading within the acceptable range.
- A. Turn the Control Switch to "X100" and connect a variable frequency pulse generator to test point 2 (pin 5 of the 24 pin connector).
 - B. Adjust the frequency to the value required to calibrate the unit at 65% of full scale in the meter (see Table T-1).
 - C. Adjust R28 (50 Kohm X100 calibration pot) until the meter reads 65% of full scale.
 - D. Readjust the frequency to the value required to calibrate the unit at 35% of full scale (see Table T-1) and note the meter reading.
 - E. If the reading is within the acceptable range listed in Table T-1, note the meter readings from steps C and D on a Certificate of Calibration. If not, readjust the frequency to the 65% value and adjust R28 to produce a different
 - G. Repeat steps D and E until the meter readings obtained are within the acceptable ranges for both the 65% and 35% of full scale calibration points.
 - H. Repeat steps B, C, D, E and F (if necessary) for the X10 range (using R26, the 500 Kohm X10 calibration pot), the X1 range (using R24, the 5 megohm calibration pot), and the X0.1 range (using R22, the 2 megohm X0.1 calibration pot).
 - I. If specified to be supplied with an SWGM or EWGM probe, perform an isotopic calibration:
 - 1) Connect the GM probe to the probe connector via a cable.
 - 2) Place the probe in a known ^{137}Cs radiation field in a fixed geometry.

Table T-2

Range	Field Strength (mR/h)	Acceptable Meter Reading (mR/h)	
X1000 (65%)	1300	1170.	- 1430
X1000 (35%)	700	630.	- 770
X100 (65%)	130	117.	- 143
X100 (35%)	70	63.	- 77
X10 (65%)	13	11.7	- 14.3
X10 (35%)	7	6.3	- 7.7
X1 (65%)	1.3	1.17	- 1.43
X1 (35%)	0.7	0.63	- 0.77
X0.1 (65%)	0.13	0.117	- 0.143
X0.1 (35%)	0.07	0.063	- 0.077

QC Acceptance Procedure

2.0 Calibration (cont'd)

- 3) Calibrate each range (except X1000) at 65% and 35% of full scale using the values in Table T-2 and the same calibration pots as those listed for a CPM calibration.
 - 4) Note the meter readings from Step 3 on a Certificate of Calibration.
10. The "X1000" range (internal detector) must always be isotopically calibrated.
- A. Position the unit so that the internal GM tube is in the proper fixed geometry within a known ^{137}Cs radiation field.
 - B. Turn the Control Switch to "X1000" and calibrate at 65% and 35% of full scale, using the values in Table T-2 and R30 (the 5 Kohm X1000 calibration pot).
 - C. Note the meter readings from Step 2 on a Certificate of Calibration.
11. Mark all trimpot bodies in such a way as to show the general position of the adjustment screws after calibration.

3.0 Response Control Test

1. Connect a variable frequency pulse generator to test point 2 (pin 5 of 24 pin connector) and turn the Response Switch to "FAST" and the Control Switch to "X1000".
2. Vary the frequency of the pulse generator and switch between "FAST" and "SLOW", checking the operation.
3. Do the same at "X100", "X10", "X1" and X0.1"
4. Leave the pulse generator connected.

4.0 Audio Switch Test

1. Turn the Audio Switch to "ON" and check to see that the audio functions properly on all five ranges, using the pulse generator. Also, check that the audio remains silent when the Selector Switch is on "BAT." and "HV".
2. Turn the Audio Switch to "OFF" and the Alarm On/off Switch mounted on the main PC board to "ON". Drive the meter well beyond full scale with the pulse generator and observe that a continuous tone is heard. Do this on all five ranges.
3. Disconnect the pulse generator.

5.0 Anti-saturation Circuit Test

1. Turn the Control Switch to "X100" and connect a 1000 megohm resistor across the probe connector (from +900V to ground). The meter should peg beyond full scale.
2. Replace the 1000 megohm resistor with a 2000 megohm resistor. The meter should remain at zero.
3. Remove the 2000 megohm resistor.

6.0 CPM Calibration Test

1. Leave the Control Switch on "X100" and turn the Audio Switch to "ON".
2. Connect a suitable GM probe to the instrument.
3. Place a small check source near the probe and switch to "X100", "X10", "X1", and "X0.1" in turn. You should obtain a meter reading for each range.

7.0 Wrap-up

1. Remove all test equipment from the unit. Turn the Control Switch to "OFF" and install a new +9V alkaline battery (MN-1604 or equivalent) in the appropriate battery clip (either clip if the unit has the parallel-wire battery option).
2. Complete, date, and sign a Certificate of Calibration.

QC Acceptance Procedure

BLANK PAGE

User's Manual

Appendix B

Spare Parts List

Schematic Symbol	Description	Part No
	Main PC Board Assembly	1016017
C1, C11, C26	Capacitor, 0.1 μ F, Film	9211041
C2, C22, C23	Capacitor, 0.001 μ F, 1kV, cer.	9201021
C3	Capacitor, 0.001 μ F, Film	9211021
C4, C24	Capacitor, 33 μ F, 10V tan	9233361
C5, C7, C27	Capacitor, 0.01 μ F, Film	9211031
C6	Capacitor, 0.0033 μ F, Film	9213321
C8	Capacitor, 200 pF, N750 cer.	9202011
C9, C21	Capacitor, 0.22 μ F, Film	9212241
C10, C17, C25, C28	Capacitor, 0.047 μ F, Film	9214731
C12	Capacitor, 0.47 μ F, Film	9214741
C13, C15, C18	Capacitor, 1 μ F, 35 V tan	9231051
C14	Capacitor, 2.2 μ F, 16 V tan.	9232251
C19, C20	Capacitor, 0.01 μ F, 1kV cer.	9201031
D1-D6, D9-D17, D22, D23	Diode, 1N4148	9600004
D7, D8	Rectifier, 2 kV PIV	9600001
Q1, Q3	Transistor, 2N4126	9610002
Q2	Transistor, 2N4124	9610001
Q4	Transistor, 2N5210	9610005
R1, R6, R15, R17	Resistor, 10 k, 1/4 W, 1%	8510024
R2	Resistor, 200 ohm, 1/4 W, 5%	8120004
R3	Resistor, 27 k, 1/4 W, 5%	8127024
R4	Resistor, 1.3 k, 1/4 W, 5%	8113014
R5, R7, R9, R61	Resistor, 4.7 k, 1/4 W, 5%	8147014
R8, R18, R41, R51	Resistor, 470 k, 1/4 W, 5%	8147034
R10, R53	Resistor, 200 k, 1/4 W, 1%	8520034
R11, R33	Resistor, 82.5 k, 1/4 W, 1%	8582524
R12, R13, R16, R20	Resistor, 100 k, 1/4 W, 5%	8110034
R14, R57	Resistor, 274 k, 1/4 W, 1%	8527434
R19, R37, R40, R42, R43, R44, R45, R54, R58, R59,	Resistor, 1 meg., 1/4 W, 5%	8110044
R21, R47	Resistor, 499 k, 1/4W, 1%	8549934
R22	Trimpot, 2 meg.,	9392051
R24	Trimpot, 5 meg.,	9395051
R26	Trimpot, 500 k,	9395041

1016-0-U-0896

User's Manual

Appendix B (cont'd)


Spare Parts List (cont'd)

Schematic Symbol	Description	Part No
	Battery PC Board Assembly	9420001
	Case Top Assembly	1016147
	Handle	9710002
	Meter	9400016
	Meter Window	9400011
	Meter Support Bracket	9850002
	MHV Connector	9782001
	Probe Clip	9460004
	Case Bottom Assembly	1016050
	Miscellaneous	
BT1	Battery, 9V alkaline, MN1604	9750001
	Cable, Probe, MHV-MHV, 36-inch	9801001
	Knob, function	9770003
	Knob, round w/pointer	9770001

Appendix C

The drawings listed below follow this page.

Drawing Name	Number	Rev.
Schematic Circuit Diagram, Surveyor 2000	B1016927	C
Main Board Component Location, Surveyor 2000	B9700335	
Switch Board Component Location, surveyor 2000	B9700336	

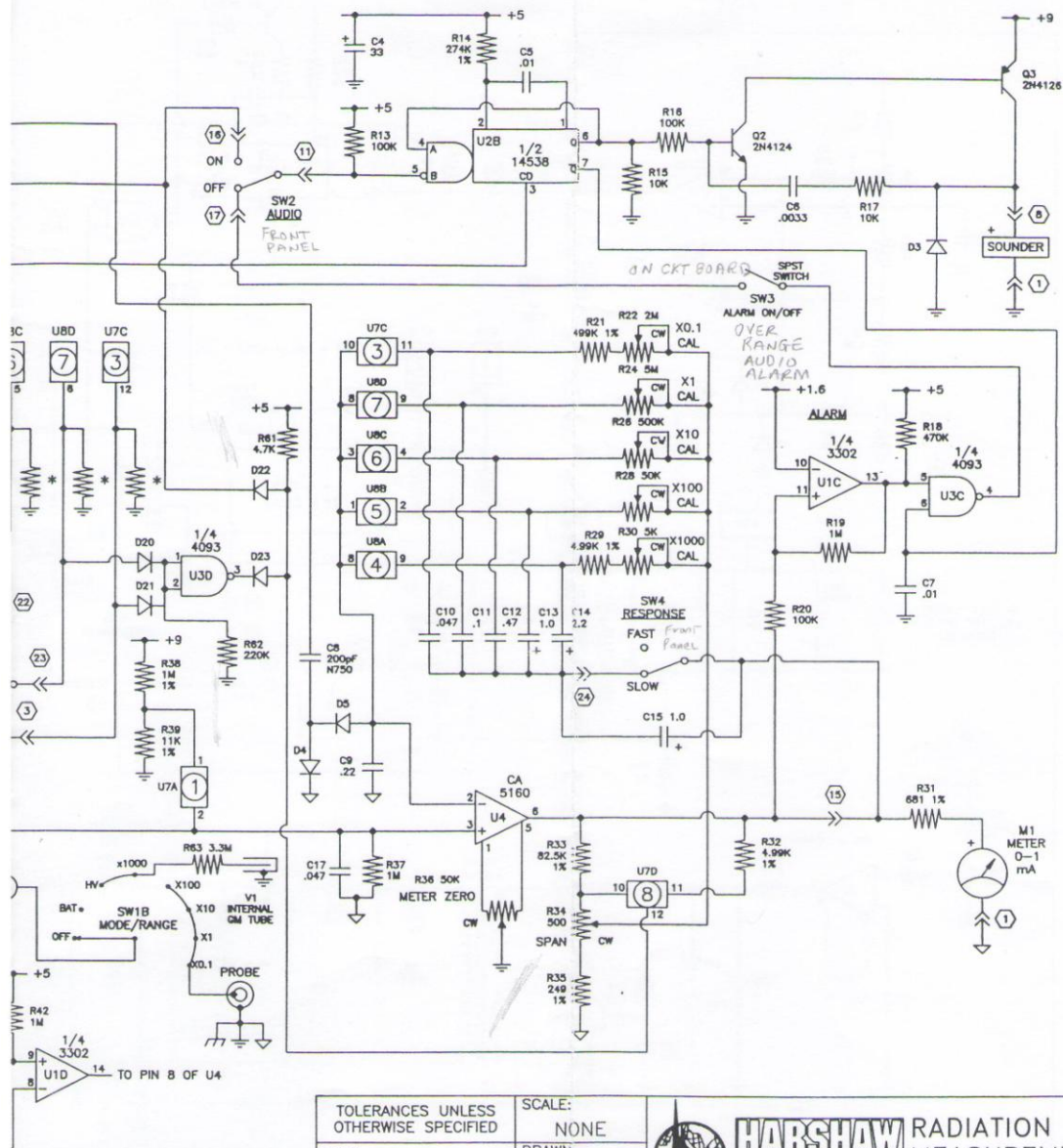

TOLERANCES UNLESS
OTHERWISE SPECIFIED

FRAC.:

.X:

.XX:

.XXX:

ANGLES:

MICRO FIN.:

DE-BURR AND BREAK
ALL EDGES

SCALE:

NONE

DRAWN

JL

DATE:

8-9-84

CHECKED:

JR

DATE:

8/3/96

DO NOT
SCALE: PRINT




**HARSHAW
BICRON** RADIATION
MEASUREMENT
PRODUCTS

SOLON, OHIO U.S.A.

Schematic Circuit Diagram

SURVEYOR 2000

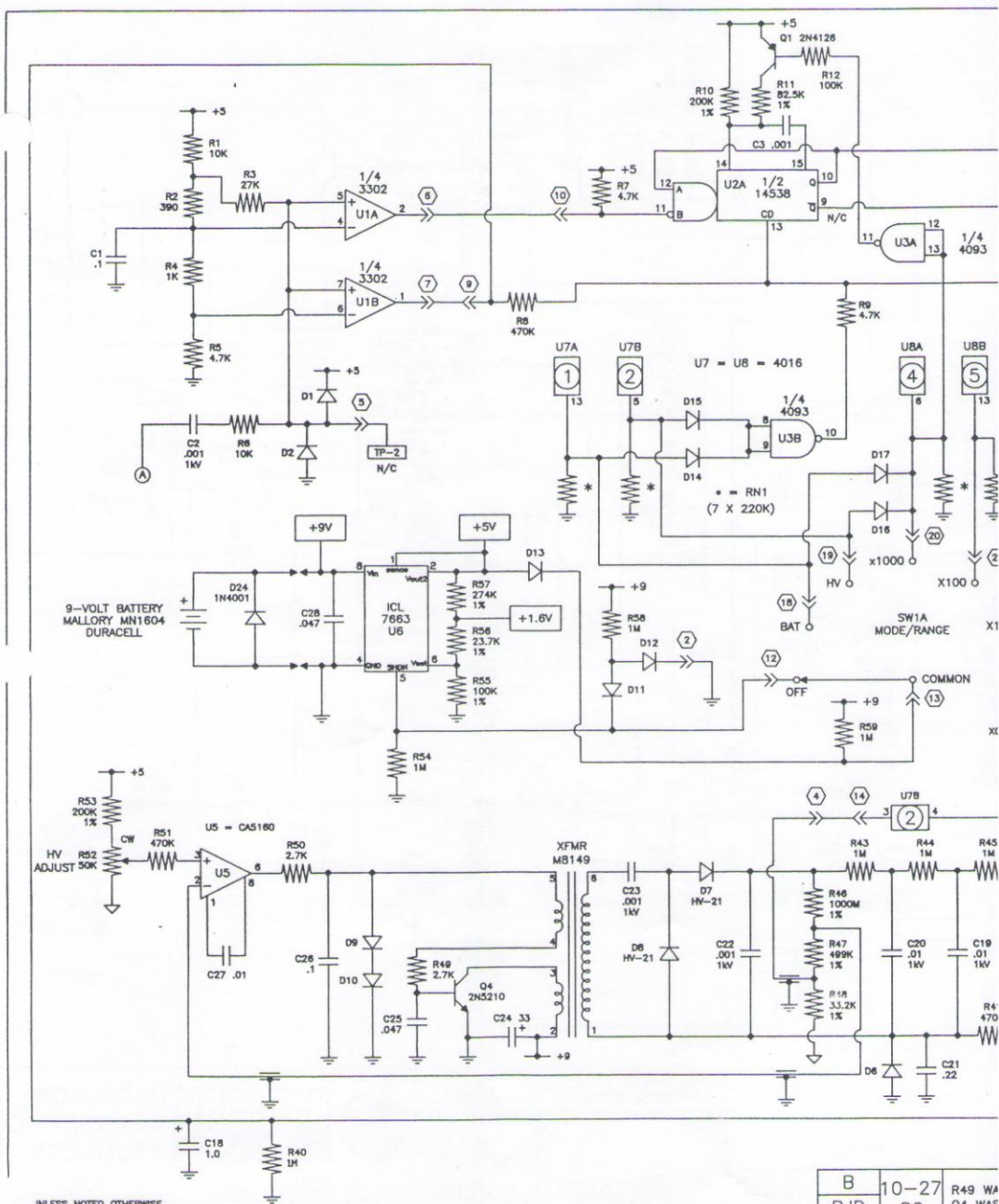

UNLESS NOTED OTHERWISE


1. ALL RESISTORS IN OHMS
2. ALL CAPACITORS IN MICROFARADS


- CONNECTOR PIN BETWEEN MAIN AND SWITCH PCB'S (CN1)
 ▶ CONNECTION BETWEEN MAIN AND BATTERY PCB'S


C	8-7	R49 WAS 1K CONVERTED TO CAD
JAR	96	
REV	DATE	DESCRIPTION
BY		

B	10-27	R49 WA Q4 WAS
RJD	86	
IN1205	A	9-4
JR	JSB	85
ECN	REV	DATE
BY	BY	

FORM

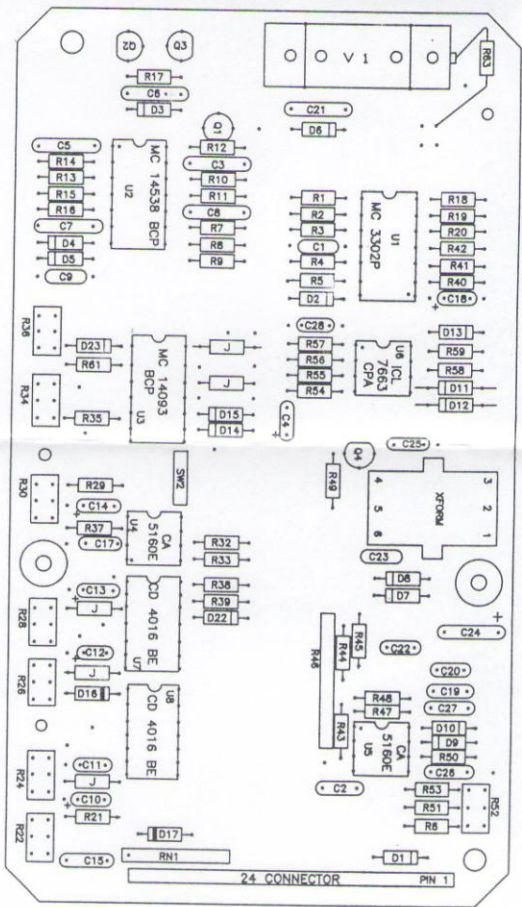

REV		DATE	DESCRIPTION	ECN	DE-BURR AND BREAK	SCALE	 HARSHAW RADIATION BICKON MEASURE PRODUCT SOLON, OHIO U.S.A.	
BY					ALL EDGES	PRINT	Main Board Component Location Surveyor 2000 9700335	
TOLERANCES UNLESS OTHERWISE SPECIFIED FRAC.: FULL .X: JAR .XX: 8-8-96 .XXX: CHECKED: JAR ANGLES: DATE: 8/3/96 MICRO FIN.: V DO NOT SCALE PRINT					DWG. SIZE: B BICK PART NUMBER: 9700335			

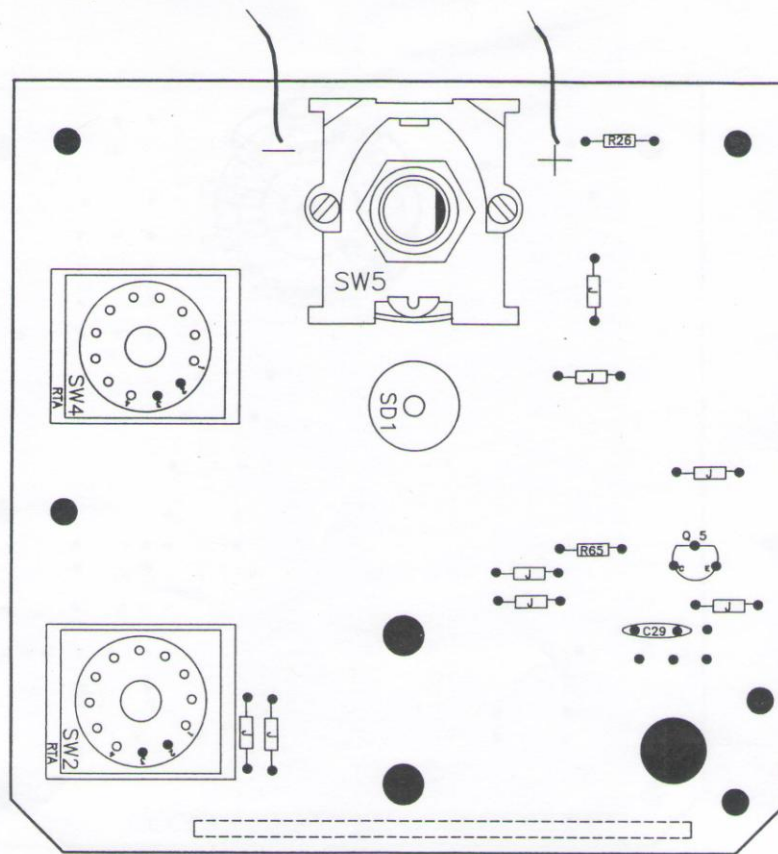

COMPONENT SIDE VIEW

REV	DATE
BY	